


Set up shading correction for color imaging in Elements

1. Go to a white area on slide (no tissue)
2. Make sure pillar (top part of microscope with bright field light) is forward and the light is on
3. Go to Shading correction tab **(A)**
4. Make sure “per optical configuration” is selected **(B)**
5. Click “Capture” **(C)**
6. In window that opens make sure “Light background (bright field)” is selected
7. Click the “Next button” approximately 8-10 times
8. Click “Finish”




Note: Only need to do this once at the beginning of imaging session!

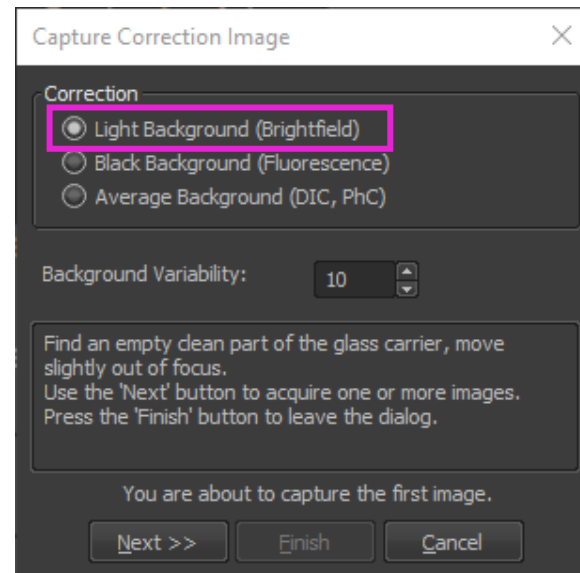


How to test shading correction

Turn on shading correction with button left of Play button

Off:

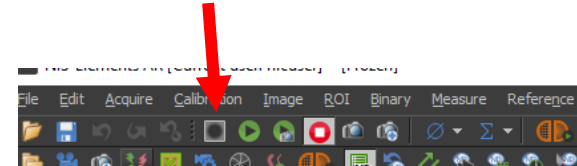

On:

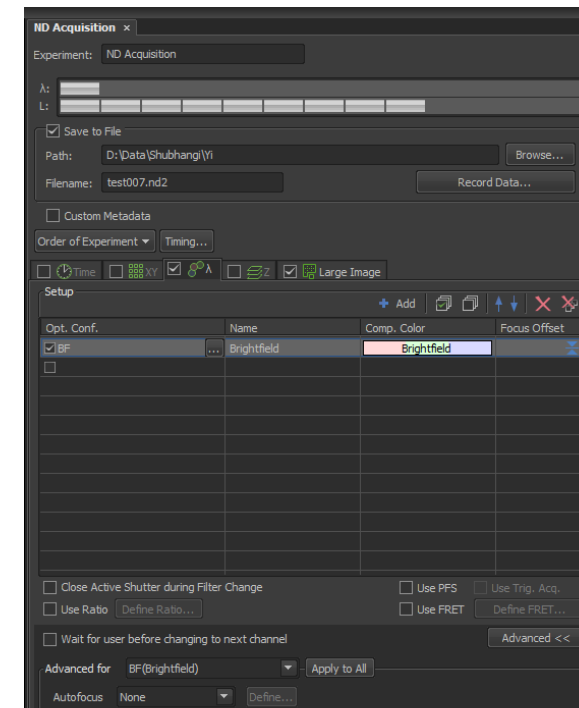

Can use ND Acquisition window to take a quick 3 x 3 image to test shading correction

- Under “λ” tab add “BF” channel
- Under “large image” set to 3x3
- Click “Run Now” on bottom right of window
- If you see the lines between the images, redo shading correction